Self-Collection with the BD Onclarity™ HPV Assay

IVD Rx Only

The FLOQSwabs® for HPV Self Collection is to be used for the collection of vaginal specimens for use with the BD Onclarity™ HPV Assay. Users should read and carefully follow instructions. If you have any questions or difficulties with the collection process described below, please contact your healthcare provider for more information. Should unexpected events not described in these instructions occur during self-collection, inform your healthcare provider.

WARNINGS: Inform your healthcare provider if you have used any substances such as vaginal creams, gels, lubricants, etc. within 24 hours prior to collection. If swab is dropped/touched in the area below the red mark, discard and request a new swab. The swab should not touch any other surface prior to inserting into vagina.


DO NOT USE: • During pregnancy or within three months after giving birth • If recently undergone gynecological surgery • In case of unusual bleeding or pelvic pain

INSTRUCTIONS:

1

Wash your hands with soap and water.

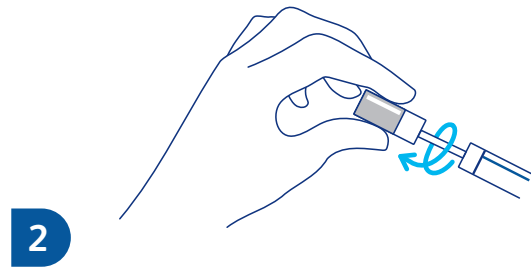2

Twist cap to break seal. Carefully pull cap off the tube to remove the swab.


Do not touch the soft swab tip or lay the swab down. If swab is dropped or touched in the area below the **red mark**, discard and request a new swab. The swab should not touch any other surface prior to inserting into the vagina.


3

Hold swab at the **red mark** with your dominant hand (the hand that you write with).

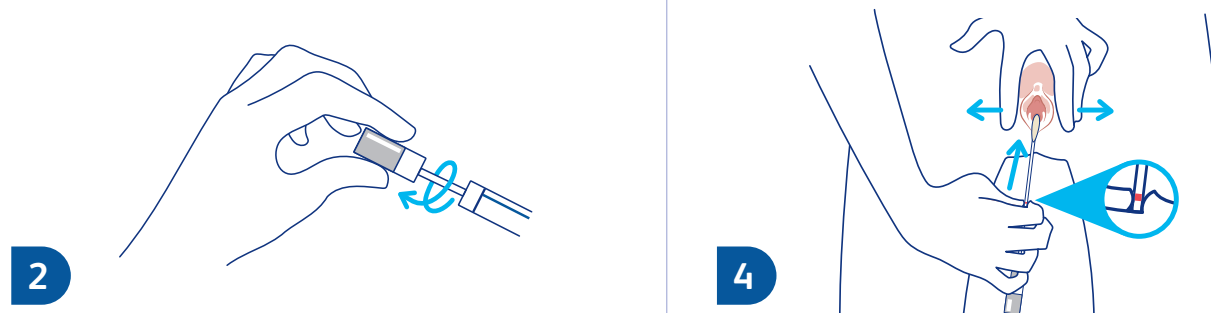
4

Standing in a comfortable position, use the free hand to spread the skin of the vaginal opening while carefully inserting the swab up to the **red mark**.

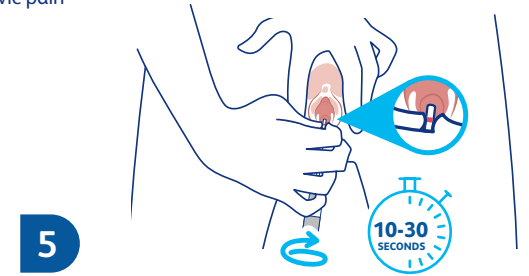
5

Gently rotate swab for 10–30 seconds and carefully withdraw the swab.

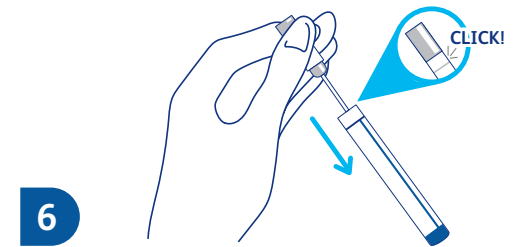6

Place swab back into the tube and securely recap.

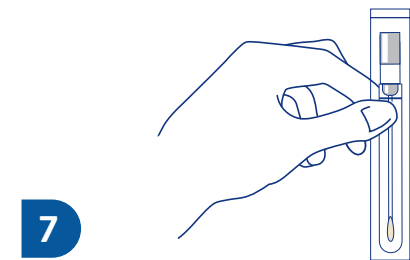7

After collection, return the recapped swab to the healthcare provider.

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
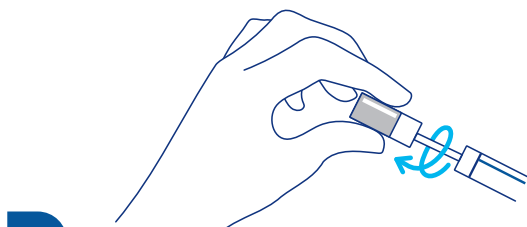
Becton, Dickinson and Company, 7 Loveton Circle, Sparks, Maryland 21152 USA

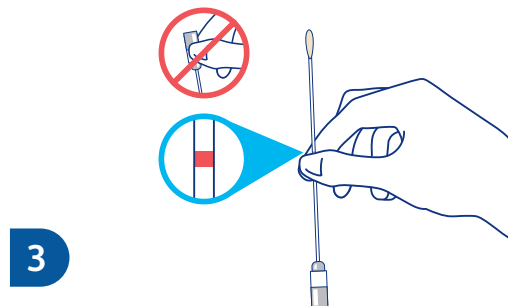BD, the BD Logo and BD Onclarity are trademarks of Becton, Dickinson and Company or its affiliates. All other trademarks are the property of their respective owners. © 2024 BD. All rights reserved.

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